




**Klausz Dental
Laboratories Ltd.**

Working **harder and smarter** for your practice!


www.klauszdentallab.com • klauszdl@uniserve.com • Tel: 416-630-6262 • Fax: 416-630-0163
40 Kodiak Crescent • Units 11-12 • Toronto, ON • Canada M3J 3G5

Mandibular Day Orthotics

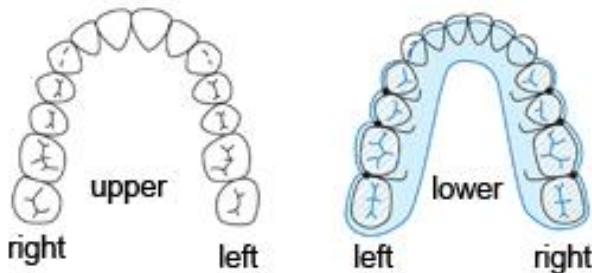


Mandibular Positioning Appliance

- * Custom PMT Base or Flexible Thermoplastic Base
- * Optional ball clasp 3-4 per side
- * Full coverage over Md posterior teeth
- * Optional anterior wrap-over incisal coverage
- * Lingual anatomical ridge 30% of Mx post cusp tips index
- * No contact of Mx buccal cusps
- * No protrusive and non-working interferences

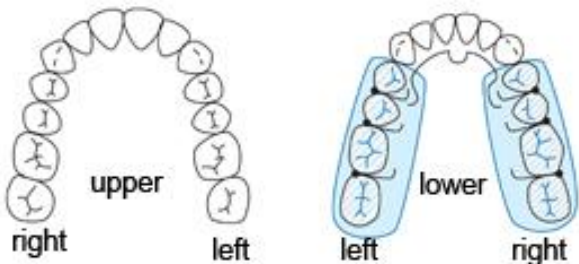


Compact Mandibular Positioning Appliance

- * Custom PMT Base
- * .032 lingual connector wire to replace major connector
- * Ball clasp 3-4 per side
- * Full coverage over Md posterior teeth
- * Lingual anatomical ridge 30% of Mx post cusp tips index
- * No contact of Mx buccal cusps
- * No protrusive and non-working interferences

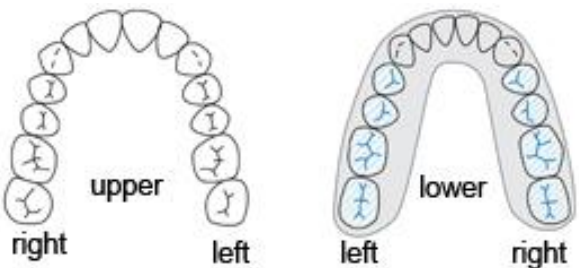

Mandibular Positioning Appliance (No Ball Claps)

- * Custom PMT Base
- * Full coverage over Md posterior teeth
- * Optional anterior wrap-over incisal coverage
- * Lingual anatomical ridge 30% of Mx post cusp tips index
- * No contact of Mx buccal cusps
- * No protrusive and non-working interferences

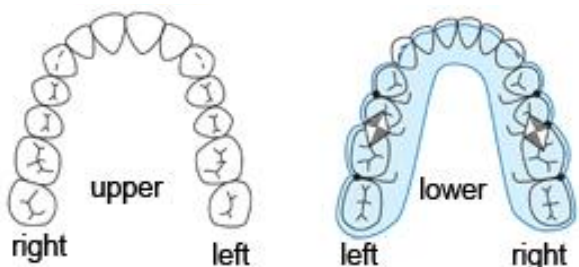



Pivot Appliance

- * Custom PMT Base
- * .032 Ball clasp 3 per side
- * Full coverage over Md posterior teeth
- * Optional anterior wrap-over incisal coverage
- * Pyramid shaped build-up over mesial buccal cusp of lower first molars and occludes against mesiolingual cusp of upper first molars.

Whatever your requirements, we can custom fabricate an orthotic to meet your specific needs!





Maxillary Night Orthotics

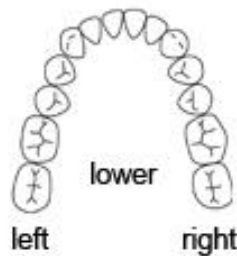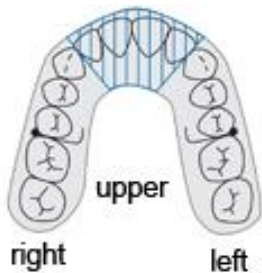



Anterior Deprogrammer

- * Custom PMT Base or Flexible Thermoplastic Base
- * Optional with ball clasp between 5-6's
- * Anterior acrylic pad contacting lower 4 incisors and extends labially and lingually to patients range of motion
- * No posterior tooth contact.

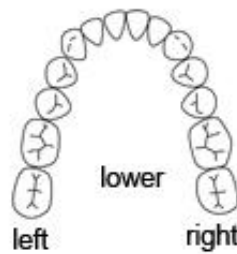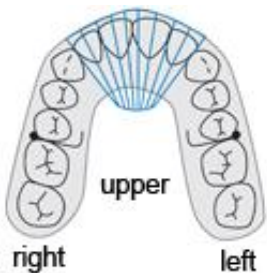



Anterior Deprogrammer with Lingual Ramp

- * Custom PMT Base or Flexible Thermoplastic Base
- * Optional with ball clasp between 5-6's
- * Anterior acrylic pad contacting lower 4 incisors and extends labially and lingually to patients range of motion
- * Solid anterior ramp contoured to lingual of lower cuspid to cuspid
- * No posterior tooth contact.

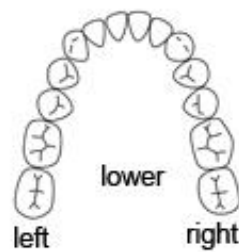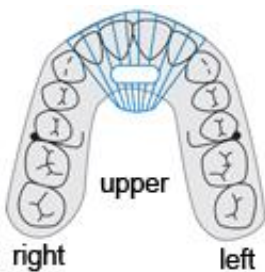

Anterior Deprogrammer with Open Lingual Ramp

- * Custom PMT Base or Flexible Thermoplastic Base
- * Optional with ball clasp between 5-6's
- * Anterior acrylic pad contacting lower 4 incisors and extends labially and lingually to patients range of motion
- * Anterior ramp with airway cut out and 2mm frame work contoured to lingual of lower cuspid to cuspid
- * No posterior tooth contact

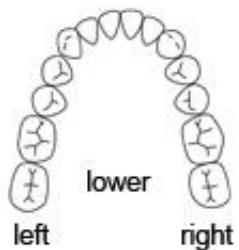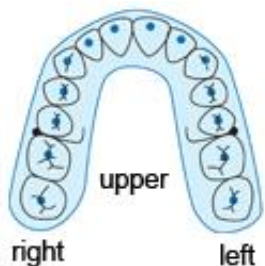

Decompressor (Flat Plane Splint)

- * Custom PMT Base or Flexible Thermoplastic Base
- * Optional with ball clasp between 5-6's
- * Mx full acrylic coverage and occlusion
- * Minimum 2mm thickness

Whatever your requirements, we can custom fabricate an orthotic to meet your specific needs!