



**Klausz Dental  
Laboratories Ltd.**

Working **harder and smarter** for your practice!

www.klauszdentallab.com • klauszdl@uniserve.com • Tel: 416-630-6262 • 1-800-410-4392  
40 Kodiak Crescent • Units 11-12 • Toronto, ON • Canada M3J 3G5

## An O.I.C.S. Specialty Laboratory *Occlusion - Implants - Cosmetics - Sleep*


# Novaloc Matrix System for Locator Abutments

The Novaloc Matrix System, with its newly developed technology, is a ready-made connective element to affix removable dentures to Locator Abutments




The snap fit connection of the Novaloc Matrix to the Locator Abutment is based on the locking mechanism of the retention inserts. These four, colour coded inserts span the functional range of the Locator Abutment and are more retentive than traditional retentive male inserts. In particular, thanks to the construction of the retention insert and the use of PEEK as a material, the Novaloc Matrix is also able to cater to situations where several extremely divergent abutment positions exist without wear and tear to the retention inserts.

- ✓ Innovative Swiss Made technology
- ✓ 4 Retentive inserts, logically colour coded according to retention values
- ✓ Retention inserts have a long lifespan and better hygiene
- ✓ Well thought-out, user friendly Novaloc tools

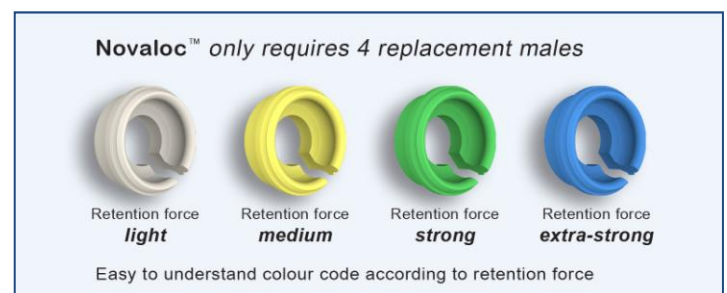

The Novaloc Matrix System offers the possibility of retention of a dental prosthesis for non-parallel implants with angles of up to 20 degrees. This means that divergences between two implants of a maximum of 40 degrees can be corrected. For this purpose, the four color-coded Novaloc retention inserts, classified by retention weight, can be used without the need for angled abutments.

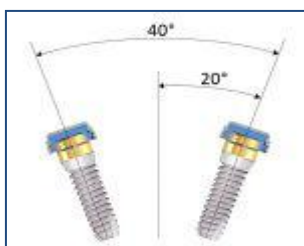

Always looking for ways to improve the quality of our work is just another way we are, "Working Harder and Smarter for your Practice!"

